ozū

MANUAL DE INSTRUCCIONES

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 www.ozulife.com

User Manual

Please read and make sure you understand the manual before use.

In particular, make sure that you fully understand the precautions.

Once the newly purchased this unit, please charge it firstly, which can protect the life of the battery. Use a 5V, 1A (or 2A) cell phone charger to charge it. When charging, the power indicator will flash, this indicates the machine is charging. After the charging is completed, the machine power indicator will be turned off.

Suitable group

People who use computers for long time

People who have poor palm blood circulation.


Women who want to maintain hand beauty.

Product features




- Ergonomic design conforms to acupuncture points and curves of hands
- Pneumatic segmented extrusion massage Full palm acupuncture massage

- Full palm acupuncture massage
- With Heat function for cold circumstance(Optional function).
- Microseismic massage increase the massage interest(Optional function).


Control Panel

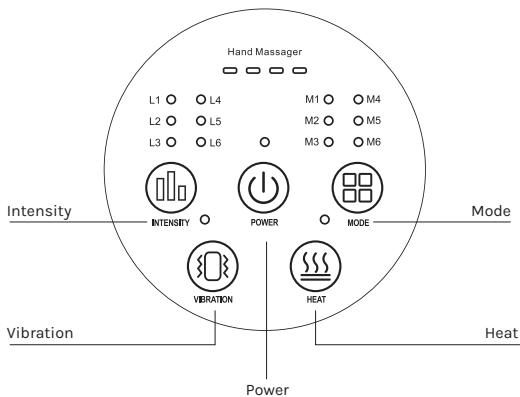Features

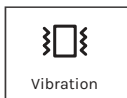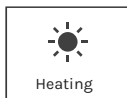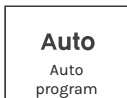Componentes

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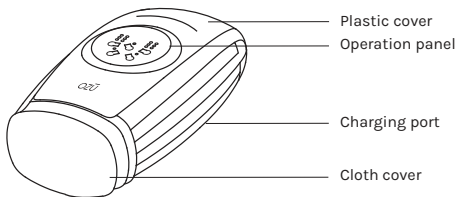Power button

Turn power on or off. Press for 2 seconds, turn on the power supply, and then press the model button, the machine will start to work. The default massage run time is 15 minutes, and the power will be turned off automatically after 15 minutes.

Intensity button

The default value is L3, with a total of 6 levels of intensity. You can select it from L3 - L4 - L5 - L6 - L1 - L2 in sequence. L1 is weakest and L6 is strongest.

Mode key

The default model is P1, with a total of 6 modes available. You can choose your favorite mode to massage, which is from P1 - P2 - P3 - P4 - P5 - P6 - P1 in sequence.

Vibration button

Press it to turn on/turn off the vibration. There are 2 modes of vibration. You can choose according to your requirement.

Heating key

Press it to turn on/turn off the heating. 5 minutes later, the inner temperature will go up obviously.

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CONDITIONS

In order to make this warranty effective, no greater requirements may be required than the presentation of the product, accompanied by the corresponding policy, duly stamped by the establishment that sold it, or the invoice, or receipt or proof, containing the specific data of the product object of the sale, depending on the case, contact the OZU offices that correspond to your country.

The consumer is responsible for the use of the product, which is not covered by the guarantee. The products contain no consumer serviceable items and do not require routine maintenance service.

THIS WARRANTY IS NOT VALID IN THE FOLLOWING CASES:

- A) When the product has been used in conditions other than normal.
- B) When the product has not been operated in accordance with the instructions for use indicated in the accompanying instructions.
- C) When the product has been altered or repaired by unauthorized entities by OZU.

In case you require parts, components, consumables and accessories, you may contact the OZU offices that correspond to your country.

TYPE OF PRODUCT: Handheld Massager Made in China

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Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
Ventis Colombia SAS. CR 10 No96 25 OF 612. NIT 901375385