
OZŪ

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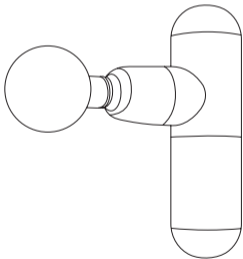COLOMBIA: Ventis Colombia SAS. CR 10 No 96 25 OF 612. NIT 901375385

INTRODUCTION

- Thank you for purchasing the OZU Portable Massage Gun.
- Please read this instruction manual carefully and keep it for future reference.
- If you require assistance, contact our support team with your model number and purchase receipt.

IMPORTANT SAFETY MEASURES

When using electrical products, especially in the presence of children, follow the basic safety precautions included in this manual:

- Children should not use this product without adult supervision.
- This device is not intended for use by individuals (including children) whose physical, sensory, or mental capabilities are reduced, or those who lack experience or knowledge, unless they have been given supervision or instructions by a person responsible for their safety.
- Children should be supervised to ensure they do not play with the device.

BEFORE USING THE PRODUCT, READ THIS MANUAL AND KEEP IT FOR FUTURE REFERENCE.

WARNINGS

To reduce the risk of burns, fires, electric shock, or injuries:

- Always unplug this product from the power source immediately after charging the battery and before cleaning it.
- Never leave the product unattended when it is plugged in. Unplug it from the power source when you have finished charging the battery and before assembling or disassembling parts.
- Do not operate it under a blanket or pillow as it can overheat.
- Do not transport this product by pulling the charging cable or use the cable as a handle.
- Never operate the product with air vents blocked. Keep air vents free of lint, hair, and similar materials.
- Do not drop or insert any objects into any openings. Do not use it outdoors.
- Do not use where aerosol (spray) products are being used or oxygen is administered.
- Ensure you unplug it from the power source when you reach for it.
- Do not unplug it with wet hands, as this could cause electric shock, electric shock-like sensations, or injuries.
- Remove the plug from the power outlet immediately after a power outage to prevent possible electric shocks, short circuits, and fires.
- Keep the cable away from hot surfaces. If the product malfunctions, stop using it immediately, and if it's charging, disconnect it. To avoid accidents or injuries, do not attempt to repair it yourself; contact the seller or manufacturer promptly.
- As the surface of this product can become hot, individuals sensitive to heat should use it with caution to prevent minor burns.
- The device's battery should only be recharged with the low-voltage safety specified on the device.

- Pay special attention to components prone to wear and please replace defective components immediately and/or keep the equipment out of use until repaired.
- Use this device only for its intended purpose as described in this manual. Cleaning and maintenance of the product should not be carried out by unsupervised children.

*If you experience discomfort during use, stop using it immediately and consult a physician. Otherwise, it may be harmful to your health.

Do not use this product if you are:

- Susceptible to electromagnetic interference, such as individuals with pacemakers or implanted medical electronic devices.
- Afflicted with malignant tumors.
- Suffering from heart diseases.
- In the early stages of pregnancy or shortly after giving birth.
- A person with an abnormal or curved spine.
- Experiencing altered perception due to a high degree of peripheral circulation failure in diabetes.
- Have scars on the skin.
- Instructed by a medical professional to rest or have health issues.
- Diagnosed with osteoporosis, severe fractures, sprains, acute muscular pain, etc.
- Running a fever of 38°C or higher.
- Experiencing acute inflammatory symptoms (fatigue, chills, changes in blood pressure, etc.) or going through an intense and debilitating stress period.
- Suffering from tendinitis.
- Being treated in medical facilities.
- Having a history of heart problems, acute illnesses, skin conditions, blood disorders such as bleeding tendencies, soft tissue diseases, heart attacks, and other physical abnormalities.



- Experiencing slow movement, physical and/or mental discomfort, or intellectual disabilities.
- Having aneurysms, bleeding issues, or taking blood thinners.
- It is not recommended to use the massager during pregnancy or menstruation. If needed, please follow your physician's recommendations.

Precautions

- This is an electronic household product for personal care and is not a substitute for medical equipment. It is not suitable for diagnosis, treatment, or commercial use.
- Before use, remove accessories from the massage area that may affect the massage.
- Stop using the product if you experience pain or discomfort in the massage area.
- If the product is not working correctly during use, turn it off immediately, disconnect it if it's charging, and contact customer support. Do not attempt to disassemble or repair the product.
- When using it in specific environments, complement its use with relevant safety regulations.
- Pay attention to electrical consumption safety. It is recommended not to plug it into the same power outlet as other devices simultaneously.
- Keep this product away from water, fire, or corrosive environments.
- Do not use this product in one place for an extended period to prevent skin damage and muscle pain due to excessive massage.
- Do not use this product while concentrating on specific operations or while driving.
- Do not use the product on areas that have recently undergone surgery, or if there are discomforts such as skin damage, redness, inflammation, congestion, etc.

- When storing the product, place it out of reach of children. It is prohibited to use it on children under 6 years of age. Children aged 6 to 12 should use the product under the supervision of a parent or the person responsible for their safety.
- Avoid using the product until one hour after your last meal.
- Do not massage directly on areas such as the head, elbows, knees, and chest.
- Do not use this product in a humid environment, and keep it away from water or other fluids.
- If the product fails, do not attempt to repair it yourself, and contact our customer service.
- It is recommended to use the product for periods not exceeding 30 minutes. Afterward, take a break.
- This device should not be used within 20 centimeters of any injury until at least 90 days have passed. Always follow your doctor's instructions.
- Do not use the device on sensitive areas such as the head, face, spine, shins, back of the feet, back of the hands, and other areas with less muscle covering the bones.

PRODUCT DIAGRAM

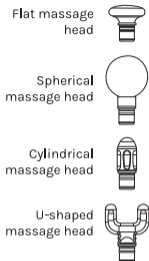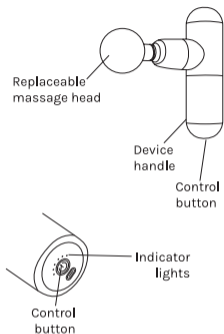OPERATING INSTRUCTIONS

Power On / Off

- Press and hold the control button for 3 seconds to turn on the device.
- Press and hold the control button for 3 seconds to turn off the device.

Massage Speed Selection

The device has indicator lights for the current motor speed. If one or more lights are on, the device is powered on.

Press the control button once more to increase the massage speed. When you reach the sixth speed, pressing the button again will return it to the first speed.

Massage Speed Levels	
One light on:	First speed (2200 rpm)
Two lights on:	Second speed (2400 rpm)
Three lights on:	Third speed (2600 rpm)
Four lights on:	Fourth speed (2800 rpm)
Five lights on:	Fifth speed (3000 rpm)
Six lights on:	Sixth speed (3200 rpm)


Timer

This device automatically shuts off after 15 minutes of use. You can turn it on again if you wish to continue using it. Remember not to use it for more than 30 minutes in a single session.


Safety Auto-Off

To protect the life of your motor, this device has an automatic safety shutdown system that activates when the head remains motionless for 5 seconds.


REPLACING THE MESSAGE HEAD

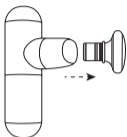Remove the message head from the device

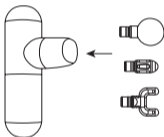
Install the message head you want to use

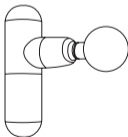Once the message head is securely in place, installation is complete

BATTERY RECHARGING

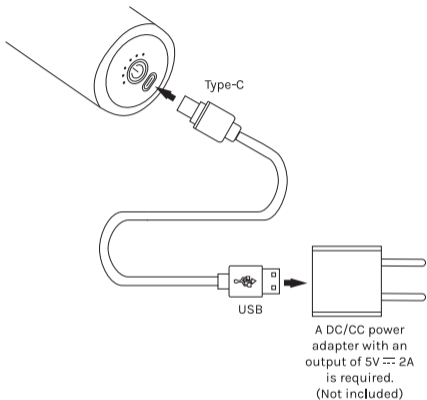- The indicator light flashes while the battery is charging and stops flashing when it is fully charged.
- To protect the battery's life, charge it as soon as the indicator light starts flashing.
- Do not use the device while the battery is charging.

PRODUCT MAINTENANCE

- Disconnect the device from the power source before cleaning it.
- Clean with a clean dry cloth.
- Do not clean the control panel with a wet cloth.
- Do not clean the product's surface with detergents or use any cleaning products.
- When the product is not in use, store it in a dry and cool place.

BOX CONTENTS

- 1 x Massage Gun
- 1 x Type-C USB Charging Cable
- 1 x Spherical Massage Head
- 1 x U-shaped Massage Head
- 1 x Cylindrical Massage Head
- 1 x Flat Massage Head
- 1 x User Manual

SPECIFICATIONS

Model:	UM-023/UM-024
Name:	Portable Massage Gun
Input:	5V $\overline{\text{---}}$ 2A
Timer:	15 min
Weight:	0.363 KG
Product type:	Massager
Country of origin:	Made in China

OZU Warranty

OZU guarantees this product for a period of 24 MONTHS from the date of purchase. Defects in materials and workmanship will be repaired or replaced by OZU, provided valid proof of purchase is presented within the specified period. The use of unauthorized replacement parts will void this warranty. The user must strictly comply with all the instructions in the instruction manual and must refrain from any action or use that is considered inappropriate or that is warned in it. If the product is used in a country where it should not be used, voltage conversion may be required. OZU is not responsible for power failures due to incorrect voltage supply to the device.

Warranty conditions

The warranty begins on the date of purchase and expires at the end of the period indicated above. A valid proof of purchase is required. OZU's warranty applies if the product has been handled properly for its intended use and used in accordance with the operating instructions at all times.

The warranty does not cover losses caused by nature, including but not limited to data loss or business loss. Some warranty exclusions apply.

Warranty Terms

The warranty begins on the day of purchase and expires at the end of the period indicated above. In the event that the merchant from which it was purchased is not an authorized OZU channel (for example, on online auction sites), the warranty period is considered to start three months after the date of manufacture indicated on the product or on the batch or serial number. Should a defect occur due to faulty workmanship or materials within the warranty period, OZU will repair it free of charge. Where repair is impossible or deemed uneconomical, OZU may agree to replace the product, at OZU's discretion.

The OZU warranty will apply as long as the product has been used correctly according to the user manual provided and the purpose it was produced for. All expenses derived from transportation to and from the service center will be borne by the customer.

Exclusions

The warranty does not cover losses due to natural causes, including but not limited to loss of data or business.

The causes of warranty cancellation are as follows:

- The purchase document was modified or altered in any way.
- The serial number, batch, or model of the product is modified, altered, removed, or made illegible.
- A repair, modification, or alteration has been carried out by unauthorized personnel or unauthorized technical service.
- The product was used for commercial purposes or other than domestic use.
- The defect was caused by abuse or misuse of the product or by environmental conditions that do not agree with what is recommended for its operation.
- The defect was caused by connection to peripheral equipment, additional equipment, or accessories other than those recommended by OZU.
- The unit has been damaged, including, but not limited to damage caused by pets, lightning, abnormal voltage, water or fire, natural disasters, or transportation accidents.
- Non-functional parts are considered consumable parts by their nature and are not included in the warranty.
- The product does not work correctly because it was not designed, manufactured, approved, and/or authorized for use in the country where it was used.
- Excesses or drops in electrical voltage that imply use in abnormal conditions, as well as defects in the electrical installation.
- Adulteration, absence, or amendment of the serial number or date of manufacture or of the label that contains the model, as appropriate.

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- Damage caused by impacts, falls, improper transport and storage, or improper treatment, as well as damage due to tremors, earthquakes, tidal waves, floods, lightning, humidity, excessive dust environments, etc.
- The operation of the product in conditions not prescribed in the user manual.
- Technical defects resulting from problems in the incorrect installation of the product.
- If applicable, the use of products, equipment, systems, spare parts, accessories, applications, external connections, installations, or connectors not supplied or recommended in the user manual.
- OZU brand products that are not imported by Ventis SpA for the Chilean market, by Ventis Perú SAC for the Peruvian market, by Ventis Colombia SAS for the Colombian market, or by Ventis Retail S.A.P.I. de C.V. for the Mexican market.

Customer Service

We recommend that you carefully read the product's operating instructions before contacting OZU or the vendor. To obtain service within the warranty period, please contact the corresponding Authorized Technical Service.

In order to assist you more efficiently, please have your purchase receipt or invoice ready when contacting OZU, indicating the date of purchase, vendor, model, and serial number, production batch, or production date as indicated on the product.

support@ozulife.com
www.ozulife.com

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CHILE: Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4