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ozū

USER MANUAL

# EYE MASK WITH HEAT AND VIBRATION MASSAGE

UM-002

WELL-BEING SERIES

## IMPORTANT SAFETY INSTRUCTIONS

Dear customer, we thank you deeply for your purchase of this device. **Please read this User Manual carefully and keep it for further reference.** When using the product, **all** instructions and basic precautions should **always** be followed, including the following:

### DANGER

#### To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

### WARNINGS:

#### To reduce the risk of burns, fire, electric shock, or injuries:

1. The product should never be left unattended when plugged in. Unplug from outlet when the device is not use and before switching in or out parts.
2. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Do not carry this product by its supply cord or use the cord as a handle.
4. Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, or similar.
5. Do not drop or insert any object into any opening.
6. Do not use outdoors.
7. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8. Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
9. Remove the power plug from the socket immediately after a power failure. Otherwise, it may cause electric shock, short circuit and fire.
10. Keep the cord away from heated surfaces. If a fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. Please do not repair it by yourself, please contact the seller or manufacturer in time.
11. As the surface of this product may get hot, people who are sensitive to heat should pay attention when using it to avoid minor burns.
12. The appliance must only be supplied at safe low voltage corresponding to the marking on the appliance.
13. The appliance is only to be used with the power supply unit provided with the appliance.
14. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
15. Special attention must be paid to components most susceptible to wear and tear. Replace defective components immediately and/or keep the equipment out of use until repair.

16. Use this appliance only for its intended use as described in this manual.
17. Children should be supervised to ensure that they do not play with the appliance.
18. Cleaning and user maintenance shall not be made by children without supervision.

**\*If you feel any abnormality during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.**

### **The following must not use this product:**

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are susceptible to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is early in pregnancy or shortly after giving birth.
6. A person who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of failure in the peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who have been medically prescribed to rest and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute headache.
11. People with 38 °C temperature (fever).
12. Acute inflammatory symptoms, (tiredness, chills, changes in blood pressure, etc.) or going through a time of intense, debilitating stress.
13. People with tendonitis.
14. Others treated at medical facilities.

### **Precautions**

1. This product is a household electronic personal care product, not a substitute for medical equipment, and not suitable for diagnosis, treatment and other commercial purposes.
2. Before use, please remove the accessories on the massage part or other items that may affect the massage.
3. Please stop using if there is pain or discomfort in the massage area when using the product.
4. When using the product, if there is any abnormality in the product, please turn off the product immediately, cut off the power supply and contact customer service, and do not disassemble and repair it by yourself.
5. When using in other specific environments, please comply with safety regulations.
6. Please use the adapter that matches the input voltage of the product to connect the product, and pay attention to the safety of the input electricity. It is recommended not to connect at the same outlet simultaneously with

other devices.

7. Please keep this product away from water, fire and corrosive environments.
8. Do not use this product in one place for a long time to prevent skin damage and muscle pain caused by excessive massage.
9. Do not use this product while focusing on precise operations or while driving.
10. Do not use it if the massage part has undergone surgery, or if there is discomfort such as skin damage, redness, inflammation, congestion, etc.
11. Patients with heart disease or medical history, acute disease, skin disease, malignant tumor, blood disease such as bleeding tendency, osteoporosis, soft tissue disease, stroke and other physical abnormalities should not use this product.
12. People with slow movement, physical and mental discomfort or intellectual disabilities should not use it. If you need to use it, please follow your doctor's advice.
13. It is not recommended to use the massager during pregnancy and menstruation. If you need to use it, please follow your doctor's advice.
14. When storing the product, please place it out of the reach of children. Children under the age of 6 are prohibited to use, and children aged 6-12 must use it under the guidance and supervision of their parents or persons with the ability to supervise.


## Important considerations

Do not use the device in the following situations:

- If you are told to stop exercising.
- If you suffer from diseases such as the following: thrombosis, severe aneurysm, acute venous aneurysm, various dermatitis, various skin infections, etc.
- If you cannot operate the mask yourself.
- If there is any discomfort during use, stop using it immediately.
- Children under 14 years of age should not use this product and should not be near them while using it.
- If you feel dizzy, vomiting and feel any discomfort during use, stop using it immediately.
- Do not abuse the massage function.
- Do not wear the mask in the bath or shower.
- Do not use under comforters or blankets.
- Do not use the device if it is hit or dropped with force, unless it is still in good condition.

Congratulations on your purchase of this OZU Eye mask with heat and vibration. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact us and send your product model number.

## Parts of the product

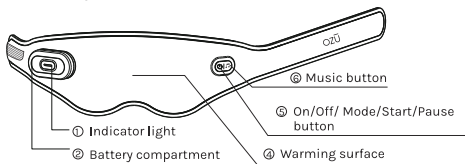- Power on/off: Press the switch button for 2 seconds to turn on and off.
- Mode change: Press the switch button for one second. There are three modes:
  - Integrated mode
  - Heating mode
  - Vibration mode

Press the switch button to change the mode in sequence. Integrated mode will start automatically when the device is running

- Music play: press the music button for 2 seconds to play music. Three different songs can be played, if you want to change music, please press the music button. Music can be paused by pressing the music button for 2 seconds while playing.
- If you need a cold compress, please put the gel ice pack in the refrigerator for 10 minutes, then take it out and put it on the eye mask.

## Warnings:

- Use only in environments where the temperature range is 5°C or 40°C.
- Pregnant women, people who have lost their sense of skin temperature, people who suffer from low-temperature burns, and people who feel unwell should consult their doctor before use. If there is any discomfort during use, please stop using it immediately.
- Each use session is limited to 15 minutes, check your eye conditions, then choose the use time you need.
- It is a normal phenomenon if your eyes turn red for a short time, but if this phenomenon continues for a long time, please consult your physician.

Indicator light	Indication	Condition
White	On and off within 5 seconds	The battery compartment is pulled into main engine and standby
White	Flash	Charging
White	Light on	Fully charged
White	Light on	Working

### Remarks:

1. The light on the power box will flash when the battery is low
2. The device can be used while it is being recharged

## Device charging

Please confirm if the battery is sufficient before using, if the battery is too low and the machine function is turned off or cannot be turned on. Please charge it before using it, be sure to unplug the power adapter after charging, to avoid charging for a long time.

1. Insert the charging plug into the charging hole on the side of the battery compartment.
2. Connect the USB socket of the power cord to a 5V= 1A charger and insert it into the power socket
3. The charging LED light flashes when charging.
4. The indicator light is continue light on when fully charged . It takes about 1 hour to fully charge it.

### NOTE:

Be sure to unplug the power cord during maintenance, but do not unplug the cord with wet hands.

## Maintenance

1. This product contains lithium battery, please do not disassemble, refit, flood the eye mask or throw it into fire, otherwise the product may be damaged.
2. If it is not used for a long time, please protect it from dust, and put the product in a dry, ventilated room, avoiding sun exposure.

## Specifications

Model No.: UM-002

Product Name: Eye Mask with Heat and Vibration Massage

Electric Specifications: 5V=1A

## Box Contents

1 x Eye mask with heat and vibration

1 x USB Type-C Cable

1 x Gel ice pack

1 x User Manual

## OZU Warranty

OZU guarantees this product for a period of 24 MONTHS from the date of purchase. Defects in materials and workmanship will be repaired or replaced by OZU, provided valid proof of purchase is presented within the specified period. The use of unauthorized replacement parts will void this warranty. The user must strictly comply with all the instructions in the instruction manual and must refrain from any action or use that is considered inappropriate or that is warned in it. If the product is used in a country where it should not be used, voltage conversion may be required. OZU is not responsible for power failures due to incorrect voltage supply to the device.

## Warranty conditions

The warranty begins on the date of purchase and expires at the end of the period indicated above. A valid proof of purchase is required. OZU's warranty applies if the product has been handled properly for its intended use and used in accordance with the operating instructions at all times.

The warranty does not cover losses caused by nature, including but not limited to data loss or business loss. Some warranty exclusions apply.

## Warranty Terms

The warranty begins on the day of purchase and expires at the end of the period indicated above. In the event that the merchant from which it was purchased is not an authorized OZU channel (for example, on online auction sites), the warranty period is considered to start three months after the date of manufacture indicated on the product or on the batch or serial number. Should a defect occur due to faulty workmanship or materials within the warranty period, OZU will repair it free of charge. Where repair is impossible or deemed uneconomical, OZU may agree to replace the product, at OZU's discretion. The OZU warranty will apply as long as the product has been used correctly according to the user manual provided and the purpose it was produced for. All expenses derived from transportation to and from the service center will be borne by the customer.

## Exclusions

The warranty does not cover losses due to natural causes, including but not limited to loss of data or business.

The causes of warranty cancellation are as follows:

- The purchase document was modified or altered in any way.
- The serial number, batch, or model of the product is modified, altered, removed, or made illegible.
- A repair, modification, or alteration has been carried out by unauthorized personnel or unauthorized technical service.
- The product was used for commercial purposes or other than domestic use.
- The defect was caused by abuse or misuse of the product or by environmental conditions that do not agree with what is recommended for its operation.
- The defect was caused by connection to peripheral equipment, additional

- equipment, or accessories other than those recommended by OZU.
- The unit has been damaged, including, but not limited to damage caused by pets, lightning, abnormal voltage, water or fire, natural disasters, or transportation accidents.
- Non-functional parts are considered consumable parts by their nature and are not included in the warranty.
- The product does not work correctly because it was not designed, manufactured, approved, and/or authorized for use in the country where it was used.
- Excesses or drops in electrical voltage that imply use in abnormal conditions, as well as defects in the electrical installation.
- Adulteration, absence, or amendment of the serial number or date of manufacture or of the label that contains the model, as appropriate.
- Damage caused by impacts, falls, improper transport and storage, or improper treatment, as well as damage due to tremors, earthquakes, tidal waves, floods, lightning, humidity, excessive dust environments, etc.
- The operation of the product in conditions not prescribed in the user manual.
- Technical defects resulting from problems in the incorrect installation of the product.
- If applicable, the use of products, equipment, systems, spare parts, accessories, applications, external connections, installations, or connectors not supplied or recommended in the user manual.
- OZU brand products that are not imported by Ventis SpA for the Chilean market, by Ventis Perú SAC for the Peruvian market, by Ventis Colombia SAS for the Colombian market, or by Ventis Retail S.A.P.I. de C.V. for the Mexican market.

## IMPORTANT

BEFORE USING THE PRODUCT, READ THE INSTRUCTIONS AND KEEP THEM FOR FUTURE REFERENCE.

## SAFETY MEASURES

When using electrical products, especially when children are present, basic safety precautions in the manual should be followed.

## Customer service

We recommend that you carefully read the product's operating instructions before contacting OZU or the vendor. To obtain service within the warranty period, please contact the corresponding Authorized Technical Service. In order to assist you more efficiently, please have your purchase receipt or invoice ready when contacting OZU, indicating the date of purchase, vendor, model, and serial number, production batch, or production date as indicated on the product.

**Chile**

+56 22 840 5567

WhatsApp: +56 9 92618728

**Peru**

0800 56185

**Mexico**

+52 55 6826 9609

WhatsApp: +52 55 1591 9179

**Colombia**

+57 320 999 6287

Telephone service hours: Monday to Friday from 9:00 a.m. to 6:00 p.m.

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CHILE: Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
COLOMBIA: Ventis Colombia SAS. CR 10 No96 25 OF 612. NIT 901375385