ozū

MANUAL DE INSTRUCCIONES

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WELL-BEING SERIES

IMPORTANT SAFETY INSTRUCTIONS

Dear customer, we thank you deeply for your purchase of this device. Please read this User Manual carefully and keep it for further reference. When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER

To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNINGS:

To reduce the risk of burns, fire, electric shock, or injuries:

1. The product should never be left unattended when plugged in. Unplug from outlet when the device is not use and before switching in or out parts.
2. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Do not carry this product by its supply cord or use the cord as a handle.
4. Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, or similar.
5. Do not drop or insert any object into any opening.
6. Do not use outdoors.
7. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8. Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
9. Remove the power plug from the socket immediately after a power failure. Otherwise, it may cause electric shock, short circuit and fire.
10. Keep the cord away from heated surfaces. If a fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. Please do not repair it by yourself, please contact the seller or manufacturer in time.
11. As the surface of this product may get hot, people who are sensitive to heat should pay attention when using it to avoid minor burns.
12. The appliance must only be supplied the safe low voltage corresponding to the marking on the appliance.
13. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
14. Special attention must be paid to components most susceptible to wear and tear. Replace defective components immediately and/or keep the equipment out of use until repair.
15. Use this appliance only for its intended use as described in this manual.
16. Children should be supervised to ensure that they do not play with the

appliance.

17. Cleaning and user maintenance shall not be made by children without supervision.
18. The appliance is only to be used with the power supply unit provided or recommended by the manufacturer.
19. This product has a heating surface. People with limited sensibility to heat should pay attention when using it.

*If you feel any abnormality during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

The following must not use this product:

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are susceptible to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is early in pregnancy or shortly after giving birth.
6. A person who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of failure in the peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who have been medically prescribed to rest and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute headache.
11. People with 38 °C temperature (fever).
12. Acute inflammatory symptoms, (tiredness, chills, changes in blood pressure, etc.) or going through a time of intense, debilitating stress.
13. People with tendonitis.
14. Others treated at medical facilities.

Precautions

1. This product is a household electronic personal care product, not a substitute for medical equipment, and not suitable for diagnosis, treatment and other commercial purposes.
2. Before use, please remove the accessories on the massage part or other items that may affect the massage.
3. Please stop using if there is pain or discomfort in the massage area when using the product.
4. When using the product, if there is any abnormality in the product, please turn off the product immediately, cut off the power supply and contact customer service, and do not disassemble and repair it by yourself.
5. When using in other specific environments, please comply with safety

regulations.

6. Please use the adapter that matches the input voltage of the product to connect the product, and pay attention to the safety of the input electricity. It is recommended not to connect at the same outlet simultaneously with other devices.
7. Please keep this product away from water, fire and corrosive environments.
8. Do not use this product in one place for a long time to prevent skin damage and muscle pain caused by excessive massage.
9. Do not use this product while focusing on precise operations or while driving.
10. Do not use it if the massage part has undergone surgery, or if there is discomfort such as skin damage, redness, inflammation, congestion, etc.
11. Patients with heart disease or medical history, acute disease, skin disease, malignant tumor, blood disease such as bleeding tendency, osteoporosis, soft tissue disease, stroke and other physical abnormalities should not use this product.
12. People with slow movement, physical and mental discomfort or intellectual disabilities should not use it. If you need to use it, please follow your doctor's advice.
13. It is not recommended to use the massager during pregnancy and menstruation. If you need to use it, please follow your doctor's advice.
14. When storing the product, please place it out of the reach of children. Children under the age of 6 are prohibited to use, and children aged 6-12 must use it under the guidance and supervision of their parents or persons with the ability to supervise.

1. Product overview

- 1.1 Common name
Transcutaneous Electrical Nerve Stimulator
- 1.2 Product model
UM-003
- 1.3 Product intended use
This Stimulator is used to promote local blood circulation, relieve pains.
- 1.4 Product performance and structure
Transcutaneous Electrical Nerve Stimulator is a product that adopts modern electronic science and technology to produce low-frequency electrical pulse in order to reach the therapeutic aim. It is to be used by adult over 18 years old. The product is mainly composed of host and conductive adhesive skin electrodes (referred to as electrode patches for short) as shown below:

The host: User can select appropriate treatment mode through function buttons on the host and know about the current status through LCD at any time. The descriptions of the function buttons are as below:

ON/OFF: Power on/ Power off

↻ : Change the mode

CH: CH1/CH2 Output Channels Selection

+ button: Intensity enhancement

- button: Intensity decrease

Electrode patches: The model is PG-2600. And the constituent materials are electrically conductive silicone, carbon film, PET and copper. The electrode patches can be used for about 80 times (2 times each day, 15 min each time), and the maximum time for each single use is 60 min. And it doesn't need to be sterilized.

Other accessories include connecting line (white buckle line, 1.2 m), batteries (AAA) and storage plate.

- 1.5 Product characteristics

- The product is small, exquisite and portable.
- Each channel is programmed with 6 modes. Applicable to different needs of people.
- Battery power display
- LCD display with a big screen makes the operation easy and clear.
- The contents displayed on LCD can guide the treatment; easy to operate.

2. Technical parameters and specification

- 2.1 Product specification

Product Name	Transcutaneous Electrical Nerve Stimulator
Model	UM-003
Power supply	d.c.3V (AAA LR03 battery *2) 60mA
Safety category	BF type
Product service life	≥ 3 years
Product dimensions	105.4*55*29.5 (mm)
Storage time	≤6 months

- 2.2 Condiciones ambientales

	Temperature	Humidity	Barometric pressure
Normal working condition	+5°C -+40°C	15%-93%RH	700hPa-1060hPa
Storage condition	-25°C-+70°C	0-93% RH	700hPa-1060hPa
Condiciones de transporte	-10°C-40°C	15% - 93% RH	700hPa-1060hPa

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Pulse frequency	1Hz—400Hz
Pulse width	20 μ s ~ 220 μ s
Direct current component	0V
Peak pulse voltage	<60V
Electric quantity of single pulse when the output amplitude is maximum	> 7 μ C
Maximum output energy of single pulse	\leq 300mJ
The effective value of maximum output amplitude	\leq 25V (50mA)
The peak value of output voltage when conducting open loop measurement	\leq 500V
The effect of output end open circuit and short circuit	Can endure the effect of output end open circuit and short circuit and the performance should not be weakened.
Adjustment of output amplitude	Continuous and uniform, minimum output should be no bigger than 2% of maximum output.
Rated load impedance	500 Ω with a permissible deviation of \pm 10%
Timing time of therapy stimulator	15 minutos con una desviación permitida de \pm 10%

3. Contraindications

- Pregnant, sensitive patient, Cardiovascular patients or other patient who accept the treatment of doctor should consult the doctor before using it.
- Do not use this stimulator on children or people with no ability to express their own consciousness.

4. Safety precautions

- The purpose of the warning signs and symbols in this user manual is to make sure that the user can use this product safely and correctly, and thus prevent damage to the user or other people.

Warning signs, symbols and their meanings are as follows:

Warning signs	Meaning
 Warning	Indicate when wrongly used, there is a possibility of personal casualty or serious injury.
 Caution	Indicate when wrongly used, there is a possibility of personal injury or objects get damaged.

Warning

- 1) Never use this stimulator along with other medical electronic stimulators, such as cardiac pacemaker, artificial heart-lung which are used to maintain life, and electrocardiograph. Otherwise it may cause danger.
- 2) If the user is using high frequency surgical equipment and this therapy stimulator at the same time, then there might be burn in the place where the electrode patch is pasted to or damage to the stimulator.
- 3) If this stimulator is used in the vicinity (1m) of short wave or microwave therapy stimulator, then the output of this stimulator might not be stable.
- 4) Do not use this stimulator near the heart, on the head, mouth cavity, private parts, parts with skin disease and so on.
- 5) Do not use this stimulator in places where there are high heat, inflammables and electromagnetic radiation.
- 6) Do not use this stimulator while bathing or sleeping.
- 7) Do not modify the stimulator, or else electric shock might occur.
- 8) Do not inhale or swallow small parts, which might lead to life-threatening condition.
- 9) Do not arbitrarily disassemble, repair and reform this stimulator, otherwise it might cause malfunction or electric shock accidents.

Caution

- 1) Do not move the stimulator while using, or when the user wishes to change treated area, first be sure to turn off the stimulator then use it, or it might cause strong stimulation.
- 2) If the user doesn't feel well due to stimulator anomaly, please stop using at once and consult a doctor.
- 3) If the stimulator is not to be used in a long time, please take out the battery to avoid leakage which may damage the stimulator.
- 4) Do not use this stimulator under circumstances that are beyond the range of application.
- 5) Do not use or contact with the product for a long time.
- 6) The patient is an intended operator.
- 7) Warning against servicing and maintenance while the me equipment is in use, functions the patient can safely use and, where applicable, which functions the patient cannot safely use; and maintenance the patient can perform.

5. Operation method

- 5.1 Preparations before use
 - Open battery cover and insert two AAA batteries into the battery holder of host according to polarity;


Attention:

 - 1) The host uses "AAALR03" battery which can be used for about 1 months (use one time each day, 15min each time).
 - 2) When the intensity of stimulation is very weak, it means the battery is low and a new battery should be replaced.
 - 3) When the product is not to be used for a long time, please take out the battery. Or else battery leakage may cause product malfunction.
 - 4) Please do not use batteries that are not of specified type.
 - 5) Please do not install the battery's positive and negative poles wrongly. See clearly the "+" "-" marks.
 - 6) Please replace batteries after the power has been shut off.
 - 7) Dispose the battery according to local pollution control regulation.
 - Insert one end of two-core connecting line to host socket and connect the other to electrode patch.

Attention: The connecting line can only be used to connect electrode patch and the host, it is prohibited to be used for other purposes. Do not put the connecting line around the neck in case of strangulation.

- 5.2 Get started

Tear off the protective films on electrode patches, then paste the electrode patches onto body parts in need of stimulation, including "shoulders, back, arms, legs, feet and joints (as shown below), not applicable to heart and upper back. Push the switch from OFF to ON to start up the device. It's recommended to use the device twice a day, 15min each time.

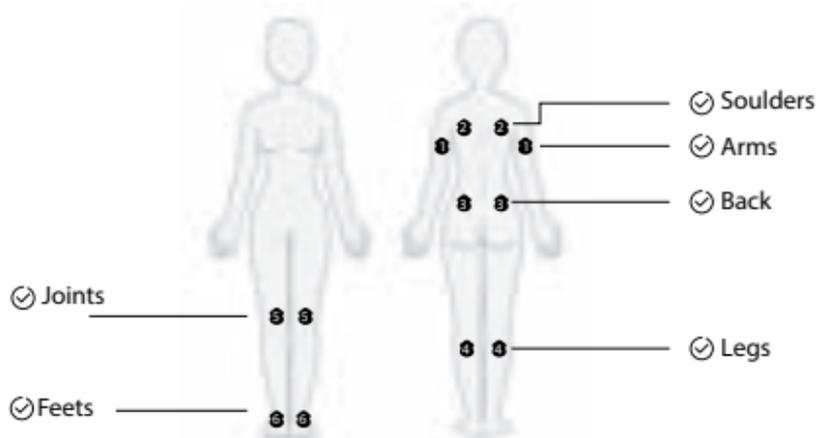Attention:

- 1) The electrode patch can be used for about 80 times (15min each time). When its viscosity can't be restored even after being cleaned for many times, please purchase new electrode patch from retailer or manufacturer.
- 2) Only use electrode patches equipped with the product. The size of electrode patch is about 55mm (width)*79mm (height)* 2.2mm (thickness). Please do not use electrode patch of other size.
- 3) Do not dispose the electrode patch casually. Please follow the local environmental requirements.
- 4) Make sure the connection between electrode patch and host is good, or else it might affect the function of the product.
- 5) One pair of electrode patches can only be used by one person.

- 5.3 Usage method of the host

1. It has two separate channels which could work individually or simultaneously. Each channel is programmed with 6 modes (12 modes in total) and 16-degree intensity.
2. Tear the protective film of electrode patches, stick them onto sore parts (Ashi points), keep relax, hold the host and turn "ON/OFF" to start.
3. After starting up, the default channel is framed "CH1" and timing 15 minutes. Press " " to change the mode while "+" "-" buttons to adjust the intensity. Press "CH" to start CH2 channel, then "CH2" on LCD screen is framed. Activate

CH2 by “ ”+“-”buttons.

4.It has safety protection function. When electrode patches are not well attached or lose effect,only 1-degree intensity is allowed. If losing effect, users are suggested to replace new patches otherwise intensity could not be strengthened. Once shift mode, intensity returns to 1. Users need to reset by the most comfortable feeling.

5.Working time of selected channel could be adjusted by long pressing “CH”, each long pressing will add 5 minutes, the longest is 30 minutes.

Each startup, the default time is 15 minutes, CH1, mode 1 and intensity “0”.

After 1 minutes without any operation, device will turn off automatically

- 5.4 After usage

After using the low-frequency therapy stimulator, please paste the electrode patches on storage plate and store in dry, ventilated place for next use.

- 5.5 Treatment drawing of common disease

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6. Maintenance and storage

- 6.1 Cleaning and maintenance

- 1)When the host is dirty, use dry soft cloth or dry towel to wipe it.
- 2)Keep the surface of electrode patches clean, avoid dust, oily matter, viscous matter, or else the viscosity of electrode patches will lower. When the surface of electrode patch is dirty, use water to clean it and dry it before use it again.
- 3)User does not need to calibrate or maintain the stimulator, and is prohibited from opening, disassembling and maintaining the stimulator arbitrarily. If there's something wrong, please contact the dealer for maintenance.
- 4)Please take out the batteries if the stimulator is not to be used for a long time.

- 6.2 Storage

- 1) Please keep the stimulator out of children's reach.
- 2) Do not store the stimulator in places where it will be exposed to direct sunlight, high temperature and moisture.
- 3) Please store the stimulator in dry and ventilated places.
- 4) Do not arbitrarily disassemble, repair and reform this stimulator, otherwise it might cause accident or malfunction.

Note: The apparatus shall be strictly disposed as per local environmental protection

7. Product disposal

In order to avoid environmental pollution, please dispose this stimulator according to local environmental requirements and do not discard casually.

8. Troubleshooting

Fault phenomenon	Reason	Solution
No display in LCD	Is the battery running out or placed inversely?	Replace new battery or put the battery correctly.
	Are there any foreign bodies in the battery holder?	Check and remove
	Is the contact between battery and battery spring piece bad?	Use suitable tool to scrape the battery spring piece clean.
The LCD displays normally, but there's no stimulation felt.	Is the two-core connecting line connected aright?	Connect the two-core connecting line correctly.
	Has the protective film on the electrode patch been torn off?	Tear off the protective film.
There is a sense of piercing pain and the skin is red.	Is the treatment intensity too high?	The suitable time is 15 minutes.
	¿La intensidad del tratamiento es demasiado alta?	Turn down the intensity properly.

<p>There is a sense of piercing pain and the skin is red.</p>	<p>Is the skin allergic to electrode patch?</p>	<p>Check whether the skin has an allergic history. If the allergy is mild, then all the user have to do is to just change the position where the electrode patches are pasted to or shorten the treatment time.If the allergy is severe, you should stop using and receive desensitization treatment first.</p>
<p>There is a sense of piercing pain and the skin is red.</p>	<p>Are the electrode patches dirty?</p>	<p>Clean the electrode patches.</p>
	<p>Are the electrode patches pasted closely to the skin?</p>	<p>Please paste the electrode patches closely to the skin.</p>
<p>The stimulation is weak.</p>	<p>Are the electrode patches pasted closely to the skin?</p>	<p>Please paste the electrode patches closely to the skin.</p>
	<p>Are the electrode patches overlapping with each other?</p>	<p>Please separate the electrode patches and then paste.</p>
	<p>Are the electrode patches dirty?</p>	<p>Clean the electrode patches.</p>


9. EMC statement

- 1.Model UM-003 need special precautions regarding EMC and need to be installed and put into service according to the EMC information provided in the accompanying document;
- 2.Portable and mobile RF communications equipment can affect models UM-003.

Warning:

- 1) The use of accessories, transducers and cables other than those specified with the exception of transducers and cables sold by the manufacturer of the models UM-003 as replacement parts for internal components, may result in increased emissions or decreased immunity of the models UM-003.
- 2) Models UM-003 should not be used adjacent to or stacked with other equipment Guidance and Manufacturer's declaration can be seen in the attachment.

9.Simbología

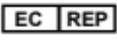Símbolo	Significado
	Batch number
	Serial number
	Date of manufacture
	Manufacturer
	Authorised representative in the European community
	Caution
	Applied part of type BF
	This product is for indoor use only
	Keep out of the sun
	Low-frequency electromagnetic radiation
	Keep away from rain

	CE mark and noticed body code
	"WEEE (Waste Electrical and Electronic Equipment)". The waste products should be handled legally.
	Product package should be recycled.
	Refer to instruction manual
<p data-bbox="134 571 217 613">IP22</p>	Dustproof waterproof level. It can prevent solid object larger than 12mm from intruding, and when tilt for 15 degrees, it can still prevent water from intruding, so no harmful effect will be created.

11. Executive standards

The product conforms to the following standards and laws:

- 1) EN 60601-1 Medical electrical equipment — Part 1: General requirements for basic safety and essential performance.
- 2) EN60601-1-11 Medical electrical equipment General requirements for basic safety and essential performance- Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.
- 3) EN60601-2-10 Medical electrical equipment - Part 2-10: Particular requirements for the safety of nerve and muscle stimulators.
- 4) EN60601-1-2 Medical electrical equipment — Part 1-2: General requirements for safety — Collateral standard: Electromagnetic compatibility — Requirements and tests.

OZU Warranty

OZU guarantees this product for a period of 24 MONTHS from the date of purchase. Defects in materials and workmanship will be repaired or replaced by OZU, provided valid proof of purchase is presented within the specified period. The use of unauthorized replacement parts will void this warranty. The user must strictly comply with all the instructions in the instruction manual and must refrain from any action or use that is considered inappropriate or that is warned in it. If the product is used in a country where it should not be used, voltage conversion may be required. OZU is not responsible for power failures due to incorrect voltage supply to the device.

Warranty conditions

The warranty begins on the date of purchase and expires at the end of the period indicated above. A valid proof of purchase is required. OZU's warranty applies if the product has been handled properly for its intended use and used in accordance with the operating instructions at all times.

The warranty does not cover losses caused by nature, including but not limited to data loss or business loss. Some warranty exclusions apply.

Warranty Terms

The warranty begins on the day of purchase and expires at the end of the period indicated above. In the event that the merchant from which it was purchased is not an authorized OZU channel (for example, on online auction sites), the warranty period is considered to start three months after the date of manufacture indicated on the product or on the batch or serial number. Should a defect occur due to faulty workmanship or materials within the warranty period, OZU will repair it free of charge. Where repair is impossible or deemed uneconomical, OZU may agree to replace the product, at OZU's discretion. The OZU warranty will apply as long as the product has been used correctly according to the user manual provided and the purpose it was produced for. All expenses derived from transportation to and from the service center will be borne by the customer.

Exclusions

The warranty does not cover losses due to natural causes, including but not limited to loss of data or business.

The causes of warranty cancellation are as follows:

- The purchase document was modified or altered in any way.
- The serial number, batch, or model of the product is modified, altered, removed, or made illegible.
- A repair, modification, or alteration has been carried out by unauthorized personnel or unauthorized technical service.
- The product was used for commercial purposes or other than domestic use.
- The defect was caused by abuse or misuse of the product or by environmental conditions that do not agree with what is recommended for its operation.
- The defect was caused by connection to peripheral equipment, additional equipment, or accessories other than those recommended by OZU.
- The unit has been damaged, including, but not limited to damage caused by

pets, lightning, abnormal voltage, water or fire, natural disasters, or transportation accidents.

- Non-functional parts are considered consumable parts by their nature and are not included in the warranty.
- The product does not work correctly because it was not designed, manufactured, approved, and/or authorized for use in the country where it was used.
- Excesses or drops in electrical voltage that imply use in abnormal conditions, as well as defects in the electrical installation.
- Adulteration, absence, or amendment of the serial number or date of manufacture or of the label that contains the model, as appropriate.
- Damage caused by impacts, falls, improper transport and storage, or improper treatment, as well as damage due to tremors, earthquakes, tidal waves, floods, lightning, humidity, excessive dust environments, etc.
- The operation of the product in conditions not prescribed in the user manual.
- Technical defects resulting from problems in the incorrect installation of the product.
- If applicable, the use of products, equipment, systems, spare parts, accessories, applications, external connections, installations, or connectors not supplied or recommended in the user manual.
- OZU brand products that are not imported by Ventis SpA for the Chilean market, by Ventis Perú SAC for the Peruvian market, by Ventis Colombia SAS for the Colombian market, or by Ventis Retail S.A.P.I. de C.V. for the Mexican market.

IMPORTANT

BEFORE USING THE PRODUCT, READ THE INSTRUCTIONS AND KEEP THEM FOR FUTURE REFERENCE.

SAFETY MEASURES

When using electrical products, especially when children are present, basic safety precautions in the manual should be followed.

Customer service

We recommend that you carefully read the product's operating instructions before contacting OZU or the vendor. To obtain service within the warranty period, please contact the corresponding Authorized Technical Service.

In order to assist you more efficiently, please have your purchase receipt or invoice ready when contacting OZU, indicating the date of purchase, vendor, model, and serial number, production batch, or production date as indicated on the product.

Chile

+56 22 840 5567

WhatsApp: +56 9 92618728

Peru
0800 56185

Mexico
+52 55 6826 9609
WhatsApp: +52 55 1591 9179

Colombia
+57 320 999 6287

Telephone service hours: Monday to Friday from 9:00 a.m. to 6:00 p.m.

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CHILE: Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
COLOMBIA: Ventis Colombia SAS. CR 10 N°96 25 OF 612. NIT 901375385