ozū

MANUAL DE INSTRUCCIONES

# MASAJEADOR SHIATSU CON CALOR

UM-008

WELL-BEING SERIES

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Ventis Colombia SAS. CR 10 No96 25 OF 612. NIT 901375385

 [www.ozulife.com](http://www.ozulife.com)

Dear customer we thank you deeply for your purchase of this device  
Please read this User manual carefully and keep it for further reference

## IMPORTANT SAFETY INSTRUCTION

When using the product, all instructions and basic precautions should always be followed, including the following:

**DANGER** - To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

## WARNINGS:

To reduce the risk of burns, fire, electric shock, or injuries:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when the device is not use and before switching in or out parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Do not carry this product by its supply cord or use cord as a handle.
- 4) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, or similar.
- 5) Do not drop or insert any object into any opening.
- 6) Do not use outdoors.
- 7) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 8) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 9) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
- 10) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 11) As the surface of this product may get hot, people who are sensitive to heat should pay attention to using it to avoid minor burns.
- 12) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 13) The appliance is only to be used with the power supply unit provided with the appliance.
- 14) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- 15) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 16) Use this appliance only for its intended use as described in this manual.
- 17) Children should be supervised to ensure that they do not play with the appliance.
- 18) Cleaning and user maintenance shall not be made by children without supervision.

\*If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.


1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.
11. People with 38 °C temperature (fever).
12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.
13. People with tendonitis.
14. Others treated at medical facilities.

## NAME AND FUNCTION OF THE MACHINE PARTS

**Front**

Massage Node

Dc connector

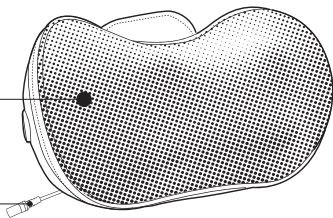**Back**


Power Adaptor

Car Adaptor
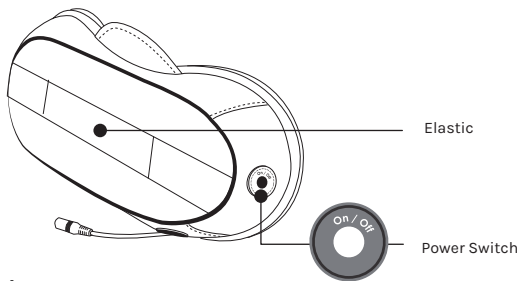
## How to use

1. Power button: Press this button, the massager will enter "starting state" and the massager will start the kneading massage function.
2. Press this button for the second time and the massage direction will switch from clockwise to counterclockwise. At the same time the heating function will activate and a red light will turn on in the massage head.
3. Press this button for the third time, the massager will shut down all massaging actions and enter the shutdown state.
4. In the startup state, the default timing time of the product is 15 minutes. The massager shuts off automatically when the timer is up and enters the shutdown state.
5. The massager is powered by power supply and can be used only when connected to the power supply.

## MASSAGE POSITION

Neck

Back

Calf

Palm

Thigh

## TO CLEAN AND MAINTAIN THE PRODUCT

### 1. Cleaning

Before you clean, please make sure that the switch is in the off location, and the power plug has been pulled out from the socket.

Please use dry cloth to clean the product.

### 2. Maintenance

Keep the product clean and far from the wet environment.

## ATTENTION

Please do not use benzene or thinner and other aggressive detergents to clean the product or it may fade.

## SPECIFICATIONS

Model: UM-008

Product Name: Shiatsu Massager with Heat

Electric Specifications: 12V $\equiv$ 2A

AC/DC Adapter

Input: 100-240V~ 50/60 Hz 0.5A

Output: 12V $\equiv$ 2A

Power consumption: 24W

Default Time: 15mins

## Box Contents

1 x Massager

1 x AC/DC Adapter

1 x Car adapter

1 x User manual

## WARRANTY

OZU warrants this product for a period of 12 months in all its parts, components and handmade labor against any manufacturing and operating defects from the date of delivery to the consumer.

OZU agrees to replace any defective part or component at no additional cost to the consumer.

## CONDITIONS

In order to make this warranty effective, no greater requirements may be required than the presentation of the product, accompanied by the corresponding policy, duly stamped by the establishment that sold it, or the invoice, or receipt or proof, containing the specific data of the product object of the sale, depending on the case, contact the OZU offices that correspond to your country.

The consumer is responsible for the use of the product, which is not covered by the guarantee. The products contain no consumer serviceable items and do not require routine maintenance service.

## THIS WARRANTY IS NOT VALID IN THE FOLLOWING CASES:

- A) When the product has been used in conditions other than normal.
- B) When the product has not been operated in accordance with the instructions for use indicated in the accompanying instructions.
- C) When the product has been altered or repaired by unauthorized entities by OZU.

In case you require parts, components, consumables and accessories, you may contact the OZU offices that correspond to your country.

## TYPE OF PRODUCT: Foot massager Made in China

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Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4