ozū

MANUAL DE INSTRUCCIONES

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UM-009

WELL-BEING SERIES

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WELL-BEING SERIES

IMPORTANT SAFETY INSTRUCTIONS

Dear customer, we thank you deeply for your purchase of this device. Please read this User Manual carefully and keep it for further reference. When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER

To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNINGS:

To reduce the risk of burns, fire, electric shock, or injuries:

1. The product should never be left unattended when plugged in. Unplug from outlet when the device is not use and before switching in or out parts.
2. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Do not carry this product by its supply cord or use the cord as a handle.
4. Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, or similar.
5. Do not drop or insert any object into any opening.
6. Do not use outdoors.
7. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8. Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
9. Remove the power plug from the socket immediately after a power failure. Otherwise, it may cause electric shock, short circuit and fire.
10. Keep the cord away from heated surfaces. If a fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. Please do not repair it by yourself, please contact the seller or manufacturer in time.
11. As the surface of this product may get hot, people who are sensitive to heat should pay attention when using it to avoid minor burns.
12. The appliance must only be supplied at safe low voltage corresponding to the marking on the appliance.
13. The appliance is only to be used with the power supply unit provided with the appliance.
14. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
15. Special attention must be paid to components most susceptible to wear and tear. Replace defective components immediately and/or keep the equipment out of use until repair.
16. Use this appliance only for its intended use as described in this manual.

17. Children should be supervised to ensure that they do not play with the appliance.
18. Cleaning and user maintenance shall not be made by children without supervision.

*If you feel any abnormality during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

The following must not use this product:

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are susceptible to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is early in pregnancy or shortly after giving birth.
6. A person who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of failure in the peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who have been medically prescribed to rest and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute headache.
11. People with 38 °C temperature (fever).
12. Acute inflammatory symptoms, (tiredness, chills, changes in blood pressure, etc.) or going through a time of intense, debilitating stress.
13. People with tendonitis.
14. Others treated at medical facilities.

Precautions

1. This product is a household electronic personal care product, not a substitute for medical equipment, and not suitable for diagnosis, treatment and other commercial purposes.
2. Before use, please remove the accessories on the massage part or other items that may affect the massage.
3. Please stop using if there is pain or discomfort in the massage area when using the product.
4. When using the product, if there is any abnormality in the product, please turn off the product immediately, cut off the power supply and contact customer service, and do not disassemble and repair it by yourself.
5. When using in other specific environments, please comply with safety regulations.
6. Please use the adapter that matches the input voltage of the product to connect the product, and pay attention to the safety of the input electricity. It is recommended not to connect at the same outlet simultaneously with other devices.
7. Please keep this product away from water, fire and corrosive environments.
8. Do not use this product in one place for a long time to prevent skin damage and muscle pain caused by excessive massage.
9. Do not use this product while focusing on precise operations or while driving.
10. Do not use it if the massage part has undergone surgery, or if there is discomfort such as skin damage, redness, inflammation, congestion, etc.

11. Patients with heart disease or medical history, acute disease, skin disease, malignant tumor, blood disease such as bleeding tendency, osteoporosis, soft tissue disease, stroke and other physical abnormalities should not use this product.
12. People with slow movement, physical and mental discomfort or intellectual disabilities should not use it. If you need to use it, please follow your doctor's advice.
13. It is not recommended to use the massager during pregnancy and menstruation. If you need to use it, please follow your doctor's advice.
14. When storing the product, please place it out of the reach of children. Children under the age of 6 are prohibited to use, and children aged 6-12 must use it under the guidance and supervision of their parents or persons with the ability to supervise.
15. The device has a heating surface, people who are not sensitive to heat need to pay attention when using it.

Essential Information

About the frequency of use

1. Use as needed, within 10 minutes to 15 minutes at a time, it is recommended to use it at appropriate intervals after using it once.
2. Excessive use may cause physical discomfort and skin injury.
3. Therefore, when using it for the first time, it is necessary to observe the physical condition and changes of the skin, and shorten the use time as much as possible.

Note before use

1. When using this product, please take a personal comfortable position and completely relax your body.
2. Before use, please remove any necklaces you are wearing to avoid discomfort during massage.
3. If you feel pain or discomfort on the skin, please stop using it and seek medical attention immediately.

About the state after use

The massage area may be red and hot for a short time, which is normal and will recover after a period of time.

Product cleaning

1. To clean the cloth or PU surface, wipe with a soft dry cloth, do not rinse or wipe with water. If you need to clean the plastic and other glossy surfaces, you can use a cloth with water or 75% alcohol to wipe the smooth surface. Do not rinse with water.
2. 75% alcohol to wipe the smooth surface. Do not rinse with water.

Others

During the normal use of the product, there will be slightly sounds such as "buzz", "pop", and "boom", which belong to the sound made by the product during normal operation.


Specifications

Product name: Travel Neck Massager
 Product Model: UM-009
 Rated voltage: 5V \approx 1A
 Battery capacity: 2500mAh
 Product weight: 670g
 Product size: 26.3*26.4*11.4 cm

Box Contents

1 x Neck Massager
 1 x USB Type-C Cable
 1 x User manual

Product Information

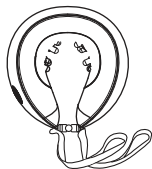Neck massager

Type-C USB
cable

Instruction Manual

Neck massage



Hands free

Heating

Product Description

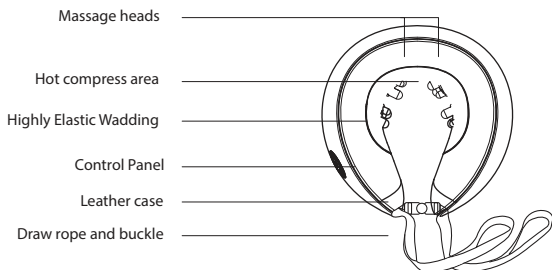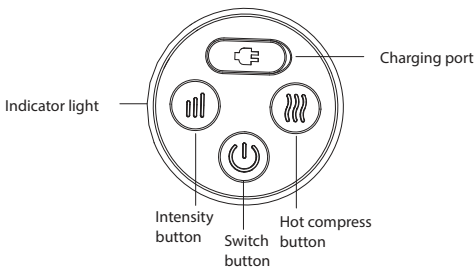Operation Instructions

Power button/Mode button 

Press and hold the center of the power button for 1.5s to turn the product on, factory mode: default second gear massage intensity + comprehensive mode, voice prompts, short press for 0.5s to adjust the massage mode (integrated mode - soothing mode - relaxation mode - kneading mode - vitality mode).

Hot compress button 

The hot compress is turned off by default. If necessary, you can short press the hot compress button to turn on the hot compress, and short press the hot compress button again to turn off the hot compress.

Strength button 

Short press the intensity button to switch the massage intensity to zero, one, two, or three gears, and cycle in ascending order;

Voice control

1. The voice mode is turned on by default. During use, you can double-click the hot compress button to turn off the voice, and double-click again to turn on the voice again.

Light color/working status	LED display	Product status
White color(when working)	Always light	Battery is normal, working normal
Red color(when working)	Always light	Low battery
White color (no charging cable plugged in when powered off)	Always light	Normal battery
Red light (no charging cable plugged in when powered off)	Always light	Low battery
Blue color(when charging)	Always light	Fully charged
Red color(when charging)	Blinking	Charging

This product cannot be used while charging. Please unplug the power adapter after charging. If the product is turned on, when the USB charging cable is plugged into the product to charge, the device will perform a shutdown operation, and a prompt sound will be emitted when the product is turned off.

Power on

Press and hold the power button for 1.5 seconds, the product is turned on, and the default mode is one (integrated mode) and the second gear is the default mode after booting.

Mode selection

Mode 1 (Comprehensive Mode)

Continuous massage for 15 minutes, forward rotation at a constant speed 10 seconds massage at rest 2 seconds and reverse rotation at a constant speed for 10 seconds cycle, the hot compress is turned off by default; if you need a hot compress, you can press the hot compress button to turn it on, and short press it again to turn off the hot compress.

Mode 2 (Soothing mode)

Continue to massage for 15 minutes, massage at a constant speed, and the hot compress is turned off by default; if you need hot compress, you can press the hot compress button to turn it on by yourself, and short press again to turn off the hot compress.

Mode 3 (Relax Mode)

Massage continuously for 15 minutes, reverse the massage at a constant speed, and the hot compress is turned off by default.

Mode 4 (kneading mode)

Reverse fast kneading for 10 seconds + slow kneading for 5 seconds cycle, the hot compress is turned off by default; if you need hot compress, you can press the hot compress button to turn it on by yourself, and short press again to turn off the hot compress.

Timer function

The default working time of each mode is 15 minutes. After 15 minutes of use, the power supply of the product will automatically cut off and shut down.

Voice broadcast

1. The voice mode is turned on by default. During use, double-click the hot compress button to turn off the voice (accompanied by a voice reminder to indicate success)

Power off

After a cycle is over/press the power button to turn off, along with the voice prompt, the state of the breathing light will be displayed during the prompt (until the voice prompt of the shutdown ends).

1. Insert one end of the Type-C charging cable into the charging port of the host.
2. Insert the other end of the Type-C charging cable into a rechargeable USB port or connect the adapter to the socket to charge the host.
 - * When the product is charging, the orange LED light will flicker during the charging process; the blue LED light will flicker when the battery is fully charged, and the light will turn off after unplugging the power cord.
 - * If you use a power adapter for power supply, you should purchase a power adapter that meets the requirements of the corresponding safety standards.

Charging time

The charging time is about 3-4 hours. When charging, the product does not support booting. When charging, the product starts a self-protection program and automatically stops working.

Finish charging

Unplug the power adapter and unplug the USB cable from the host.

Troubleshooting

If a malfunction occurs during operation, please check the following. If the fault cannot be eliminated or the product still cannot operate ideally, please cut off the power supply, and it should be repaired by professional technicians. Do not disassemble or repair this product by yourself!

Problems	Possible reason	Basic exclusion method
Can not power on	Low battery	Please charge in time
Weak massage strength	Low battery	Please charge in time
Massage area is not heated	Whether the currently selected massage mode has no heating function	Please switch the massage mode
No voice prompts during massage	Voice mode not turned on	Please follow the operation to turn on the voice prompt function

Product Disposal Requirements

Notice

- Before disposing of this product, the battery must be removed as required by regulations;
- When removing the battery, the product must be powered off;
- Batteries should be disposed of safely in accordance with laws and regulations.

Non-professionals do not disassemble by themselves!

Steps to remove the battery

1. Cut off the external power supply and remove the product shell.
2. Unplug the power supply on the battery.
3. Remove the battery.
4. Discarded batteries should be placed at a fixed recycling point.
5. Please abide by local environmental protection laws and disposal guidelines and take care of the environment.

Other Considerations

Storage

When the storage period is more than 3 months, the machine should be unpacked and recharged before use. Do not place it close to the ground, walls and roof. The room should be well ventilated to avoid strong sunlight and corrosive gas erosion.

Transportation

The product packaging is equipped with simple shockproof facilities, suitable for air, railway, road and ship transportation. Rain, snow splashes, inversions and collisions should be avoided.

Power check

Check the power supply, if the power supply voltage exceeds the specified range, do not turn it on for use.

Maintenance

Take moisture-proof, rodent-proof and pollution-proof measures. Cleaning should be done with the power off. The instrument should be placed in a safe place after use, and do not place heavy objects on the instrument. When a fault occurs, professional technicians should be invited to repair the product. Non-professionals should not disassemble the product by themselves.

Disposal and Recycling

This product should not be thrown away with household waste. When you decide to dispose of this product, please follow local environmental laws and guidelines.

OZU Warranty

OZU guarantees this product for a period of 24 MONTHS from the date of purchase. Defects in materials and workmanship will be repaired or replaced by OZU, provided valid proof of purchase is presented within the specified period. The use of unauthorized replacement parts will void this warranty. The user must strictly comply with all the instructions in the instruction manual and must refrain from any action or use that is considered inappropriate or that is warned in it. If the product is used in a country where it should not be used, voltage conversion may be required. OZU is not responsible for power failures due to incorrect voltage supply to the device.

Warranty conditions

The warranty begins on the date of purchase and expires at the end of the period indicated above. A valid proof of purchase is required. OZU's warranty applies if the product has been handled properly for its intended use and used in accordance with the operating instructions at all times.

The warranty does not cover losses caused by nature, including but not limited to data loss or business loss. Some warranty exclusions apply.

Warranty Terms

The warranty begins on the day of purchase and expires at the end of the period indicated above. In the event that the merchant from which it was purchased is not an authorized OZU channel (for example, on online auction sites), the warranty period is considered to start three months after the date of manufacture indicated on the product or on the batch or serial number. Should a defect occur due to faulty workmanship or materials within the warranty period, OZU will repair it free of charge. Where repair is impossible or deemed uneconomical, OZU may agree to replace

the product, at OZU's discretion. The OZU warranty will apply as long as the product has been used correctly according to the user manual provided and the purpose it was produced for. All expenses derived from transportation to and from the service center will be borne by the customer.

Exclusions

The warranty does not cover losses due to natural causes, including but not limited to loss of data or business.

The causes of warranty cancellation are as follows:

- The purchase document was modified or altered in any way.
- The serial number, batch, or model of the product is modified, altered, removed, or made illegible.
- A repair, modification, or alteration has been carried out by unauthorized personnel or unauthorized technical service.
- The product was used for commercial purposes or other than domestic use.
- The defect was caused by abuse or misuse of the product or by environmental conditions that do not agree with what is recommended for its operation.
- The defect was caused by connection to peripheral equipment, additional equipment, or accessories other than those recommended by OZU.
- The unit has been damaged, including, but not limited to damage caused by pets, lightning, abnormal voltage, water or fire, natural disasters, or transportation accidents.
- Non-functional parts are considered consumable parts by their nature and are not included in the warranty.
- The product does not work correctly because it was not designed, manufactured, approved, and/or authorized for use in the country where it was used.
- Excesses or drops in electrical voltage that imply use in abnormal conditions, as well as defects in the electrical installation.
- Adulteration, absence, or amendment of the serial number or date of manufacture or of the label that contains the model, as appropriate.
- Damage caused by impacts, falls, improper transport and storage, or improper treatment, as well as damage due to tremors, earthquakes, tidal waves, floods, lightning, humidity, excessive dust environments, etc.
- The operation of the product in conditions not prescribed in the user manual.
- Technical defects resulting from problems in the incorrect installation of the product.
- If applicable, the use of products, equipment, systems, spare parts, accessories, applications, external connections, installations, or connectors not supplied or recommended in the user manual.
- OZU brand products that are not imported by Ventis SpA for the Chilean market, by Ventis Perú SAC for the Peruvian market, by Ventis Colombia SAS for the Colombian market, or by Ventis Retail S.A.P.I. de C.V. for the Mexican market.

IMPORTANT

BEFORE USING THE PRODUCT, READ THE INSTRUCTIONS AND KEEP THEM FOR FUTURE REFERENCE.

SAFETY MEASURES

When using electrical products, especially when children are present, basic safety precautions in

the manual should be followed.

Customer service

We recommend that you carefully read the product's operating instructions before contacting OZU or the vendor. To obtain service within the warranty period, please contact the corresponding Authorized Technical Service.

In order to assist you more efficiently, please have your purchase receipt or invoice ready when contacting OZU, indicating the date of purchase, vendor, model, and serial number, production batch, or production date as indicated on the product.

Chile

+56 22 840 5567

WhatsApp: +56 9 92618728

Peru

0800 56185

Mexico

+52 55 6826 9609

WhatsApp: +52 55 1591 9179

Colombia

+57 320 999 6287

Telephone service hours: Monday to Friday from 9:00 a.m. to 6:00 p.m.

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CHILE: Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
COLOMBIA: Ventis Colombia SAS. CR 10 N°96 25 OF 612. NIT 901375385