ozū

MANUAL DE INSTRUCCIONES

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 www.ozulife.com

Introduction

Congratulations on acquiring your new shoulder massager UM-010. Please read this user manual before using the product and keep it for further reference.

DANGER Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

To reduce the risk of electric shock:
 Always unplug this product from the electrical outlet immediately after using and before cleaning.

To Clean and Maintain the Product

1. Cleaning

Before cleaning, please make sure that the product has been turned off and that the adaptor has been pulled out from the socket. Please use dry cloth to clean the product.

2. Maintenance

Keep the product clean and far from the wet environment.

Cómo Operar


To reduce the risk of burns, fire, electric shock, or injury to persons:

01. The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
02. Do not operate under a blanket or a pillow. Excessive heat can occur and cause fire, electric shock, or injury to persons.
03. Do not carry this product by its supply cord or use the cord as a handle.
04. Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
05. Do not drop or insert any object into any opening.
06. Do not use outdoors.
07. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
08. To disconnect, switch off all controls and remove the plug from outlet.
09. Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.
10. Keep the cord away from heated surface, otherwise it may cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
11. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
12. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
13. The appliance is only to be used with the power supply unit provided with the appliance.
14. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or persons lacking experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Product Operation Guide

1. Connect the power adaptor to the electrical outlet.
2. Press the power button to begin the massage. Heating and rolling function would also be activated.
3. Press the heating button to active heating function, or turn off the heating function.
4. Press the massage direction adjustment button to change massage direction.
5. Press the speed adjustment button to choose among 3 massage speed levels (low/medium/high).
6. Press the power button to stop the massage at any time.

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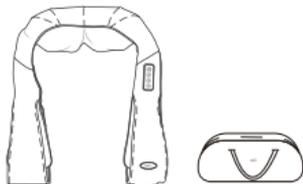Shiatsu Neck and Shoulder Massager With Carrying bag

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SPECIFICATIONS

Model: UM-010
 Product Name: Shiatsu Neck and Shoulder Massager
 Electric Specifications: 12V \approx 2A
 AC/DC Adapter
 Input: 100-240V~ 50/60Hz 0.7A
 Output: 12V \approx 2A
 Rated Power: 24W

Box Contents

1 x Massager
 1 x AC/DC Adapter
 1 x Carrying bag
 1 x User manual

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CONDITIONS

In order to make this warranty effective, no greater requirements may be required than the presentation of the product, accompanied by the corresponding policy, duly stamped by the establishment that sold it, or the invoice, or receipt or proof, containing the specific data of the product object of the sale, depending on the case, contact the OZU offices that correspond to your country.

The consumer is responsible for the use of the product, which is not covered by the guarantee. The products contain no consumer serviceable items and do not require routine maintenance service.

THIS WARRANTY IS NOT VALID IN THE FOLLOWING CASES:

- A) When the product has been used in conditions other than normal.
- B) When the product has not been operated in accordance with the instructions for use indicated in the accompanying instructions.
- C) When the product has been altered or repaired by unauthorized entities by OZU.

In case you require parts, components, consumables and accessories, you may contact the OZU offices that correspond to your country.

TYPE OF PRODUCT: Shoulder massager Made in China

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Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
Ventis Colombia SAS. CR 10 N°96 25 OF 612. NIT 901375385