ozū

MANUAL DE INSTRUCCIONES

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UM-012

WELL-BEING SERIES

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
[www.ozulife.com](http://www.ozulife.com)

## Introduction

Thank you for choosing our foot massager with heat function.

Please read the user manual carefully before use and keep it in a safe place for future reference.

If you require help, please contact our support team.
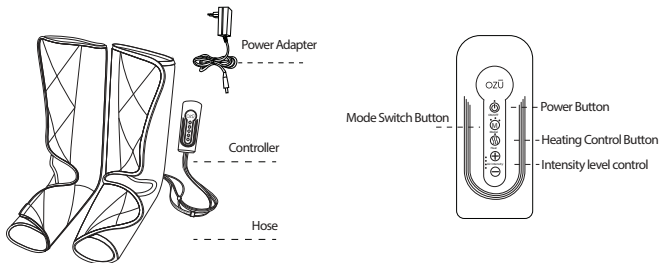
## Product Operation Guide

### How to use

Connect the adapter to the power supply and the DC power cord. Now the product is in standby mode

- **Power ON/OFF:** Press this button in the standby mode to turn on the product, and the corresponding indicator lights up. Mode1 is enabled by default, which means the heating function is automatically opened and the system starts countdown.
- **Mode button** Press this button to switch different work modes, and the corresponding function indicator is lit, There is no air compression when the indicator is off. You can use the heat function separately.

Air compression mode 1: Combination massage for feet and calves.

Air compression mode 2: Feet massage.

Air compression mode 3: Calves massage.

- **Heating control button:** The heating function is enabled by default, and the corresponding function indicator is lit, when the heating function is turned on, press this button to turn it off, and the corresponding indicator light goes out.
- **Air compression control button:** There are 4 intensity levels available for selection. Air intensity 1 is enabled by default. Short press "+" to increase the intensity level, and the corresponding intensity indicator lights up. Short press "-" button to decrease the intensity level, and the corresponding intensity indicator lights up.

## SPECIFICATIONS

Model: UM-012

Product Name: Air Compression Leg and Foot Massager

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## CONDITIONS

In order to make this warranty effective, no greater requirements may be required than the presentation of the product, accompanied by the corresponding policy, duly stamped by the establishment that sold it, or the invoice, or receipt or proof, containing the specific data of the product object of the sale, depending on the case, contact the OZU offices that correspond to your country.

The consumer is responsible for the use of the product, which is not covered by the guarantee. The products contain no consumer serviceable items and do not require routine maintenance service.

## THIS WARRANTY IS NOT VALID IN THE FOLLOWING CASES:

- A) When the product has been used in conditions other than normal.
- B) When the product has not been operated in accordance with the instructions for use indicated in the accompanying instructions.
- C) When the product has been altered or repaired by unauthorized entities by OZU.

In case you require parts, components, consumables and accessories, you may contact the OZU offices that correspond to your country.

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Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
Ventis Colombia SAS. CR 10 N°96 25 OF 612. NIT 901375385