ozū

MANUAL DE INSTRUCCIONES

# BALANZA CORPORAL INTELIGENTE

UM-019

WELL-BEING SERIES

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UM-019

WELL-BEING SERIES

Carefully read the the instructions in this user manual and reserve it for future reference.

## Warning

Please read carefully before using the scale.

Do not use on slippery floors.

Do not use when body is wet.

## Use of the scale

Your Smart Weight Scale ships with 2 AAA batteries. Install the batteries inside the battery compartment and turn on your scale. Once turned on, you will be able to measure body weight. first install the mobile applications on their mobile devices.

The smart scale can be wirelessly connected to mobile phones or tablets. For this, mobile devices need to install the mobile application.

The easiest way to download the app is to scan the QR code below, follow these instructions.

Scan the QR code and download the app "Feelfit"

After you have downloaded the app, and created an account Press the "+" on the top of the screen and select "add device" and follow the steps on the app.

Step on the scale

Pair the scale with the app



Step on the scale again and stay still until the digital number shows in a stable position

Check the IMC and weight on the app

The unit of weight can be changed from pounds (lb.) to kilograms (kg), using the mobile app. The height unit will also change accordingly (meter vs. foot/inch).

To measure the weight correctly: place the scale on a flat, hard ground, then calibrate it before measurement. The weight must show 0.0 lbs before stepping on the scale.
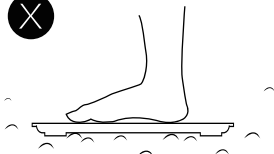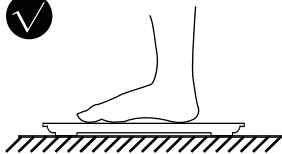
## Screen indicators

88.00  
kg

Normal Functioning

No display

Replace battery

Lo

Replace battery

EEr

Weight overload

CAL

Calibrating

## Specifications

Model: UM-019

Product Name: Smart Body Scale

Electric Specifications: 3V = (2 x 1.5V = AAA)

Batteries: 2 x 1.5V AAA

Weight range: 0.2~180kg (0.4~396lb)

Accuracy roundup: 0.1kg (0.2lb)

## Box Contents

1 x Smart Scale

1 x User manual

2 x AAA batteries

## OZU Warranty

OZU guarantees this product for a period of 24 MONTHS from the date of purchase. Defects in materials and workmanship will be repaired or replaced by OZU, provided valid proof of purchase is presented within the specified period. The use of unauthorized replacement parts will void this warranty. The user must strictly comply with all the instructions in the instruction manual and must refrain from any action or use that is considered inappropriate or that is warned in it. If the product is used in a country where it should not be used, voltage conversion may be required. OZU is not responsible for power failures due to incorrect voltage supply to the device.

## Warranty conditions

The warranty begins on the date of purchase and expires at the end of the period indicated above. A valid proof of purchase is required. OZU's warranty applies if the product has been handled properly for its intended use and used in accordance with the operating instructions at all times.

The warranty does not cover losses caused by nature, including but not limited to data loss or business loss. Some warranty exclusions apply.

## Warranty Terms

The warranty begins on the day of purchase and expires at the end of the period indicated above. In the event that the merchant from which it was purchased is not an authorized OZU channel (for example, on online auction sites), the warranty period is considered to start three months after the date of manufacture indicated on the product or on the batch or serial number. Should a defect occur due to faulty workmanship or materials within the warranty period, OZU will repair it free of charge. Where repair is impossible or deemed uneconomical, OZU may agree to replace the product, at OZU's discretion. The OZU warranty will apply as long as the product has been used correctly according to the user manual provided and the purpose it was produced for. All expenses derived from transportation to and from the service center will be borne by the customer.

## Exclusions

The warranty does not cover losses due to natural causes, including but not limited to loss of data or business.

The causes of warranty cancellation are as follows:

- The purchase document was modified or altered in any way.
- The serial number, batch, or model of the product is modified, altered, removed, or made illegible.
- A repair, modification, or alteration has been carried out by unauthorized personnel or unauthorized technical service.
- The product was used for commercial purposes or other than domestic use.
- The defect was caused by abuse or misuse of the product or by environmental conditions that do not agree with what is recommended for its operation.
- The defect was caused by connection to peripheral equipment, additional equipment, or accessories other than those recommended by OZU.
- The unit has been damaged, including, but not limited to damage caused by pets, lightning, abnormal voltage, water or fire, natural disasters, or transportation accidents.

- Non-functional parts are considered consumable parts by their nature and are not included in the warranty.
- The product does not work correctly because it was not designed, manufactured, approved, and/or authorized for use in the country where it was used.
- Excesses or drops in electrical voltage that imply use in abnormal conditions, as well as defects in the electrical installation.
- Adulteration, absence, or amendment of the serial number or date of manufacture or of the label that contains the model, as appropriate.
- Damage caused by impacts, falls, improper transport and storage, or improper treatment, as well as damage due to tremors, earthquakes, tidal waves, floods, lightning, humidity, excessive dust environments, etc.
- The operation of the product in conditions not prescribed in the user manual.
- Technical defects resulting from problems in the incorrect installation of the product.
- If applicable, the use of products, equipment, systems, spare parts, accessories, applications, external connections, installations, or connectors not supplied or recommended in the user manual.
- OZU brand products that are not imported by Ventis SpA for the Chilean market, by Ventis Perú SAC for the Peruvian market, by Ventis Colombia SAS for the Colombian market, or by Ventis Retail S.A.P.I. de C.V. for the Mexican market.

## IMPORTANT

BEFORE USING THE PRODUCT, READ THE INSTRUCTIONS AND KEEP THEM FOR FUTURE REFERENCE.

## SAFETY MEASURES

When using electrical products, especially when children are present, basic safety precautions in the manual should be followed.

## Customer service

We recommend that you carefully read the product's operating instructions before contacting OZU or the vendor. To obtain service within the warranty period, please contact the corresponding Authorized Technical Service.

In order to assist you more efficiently, please have your purchase receipt or invoice ready when contacting OZU, indicating the date of purchase, vendor, model, and serial number, production batch, or production date as indicated on the product.

### Chile

+56 22 840 5567

WhatsApp: +56 9 92618728

### Peru

0800 56185

### Mexico

+52 55 6826 9609

WhatsApp: +52 55 1591 9179

Colombia  
+57 320 999 6287

Telephone service hours: Monday to Friday from 9:00 a.m. to 6:00 p.m.

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CHILE: Ventis SpA. Camino Lo Boza 107 A09, Pudahuel, Santiago, Chile. RUT 76229287-4
COLOMBIA: Ventis Colombia SAS. CR 10 No96 25 OF 612. NIT 901375385